


**Product:** G-Rex®500-Series Gas Permeable Cell Culture Products

**REF (Catalog Numbers):**

Usage	G-Rex®500M	G-Rex®500M-CS	G-Rex®500M-TF
<b>GMP</b> Sterile Fluid Path	<b>G285500</b>	<b>G285500-CS</b>	<b>G285500-TF</b>
<b>RUO</b> Gamma Irradiated	<b>G285500-RU</b>	<b>RU05500-CS</b>	<b>RU85500-TF</b>

### 1) Description

The G-Rex®500 series products are single-use cell culture bioreactors designed for the expansion and recovery of mammalian cells.

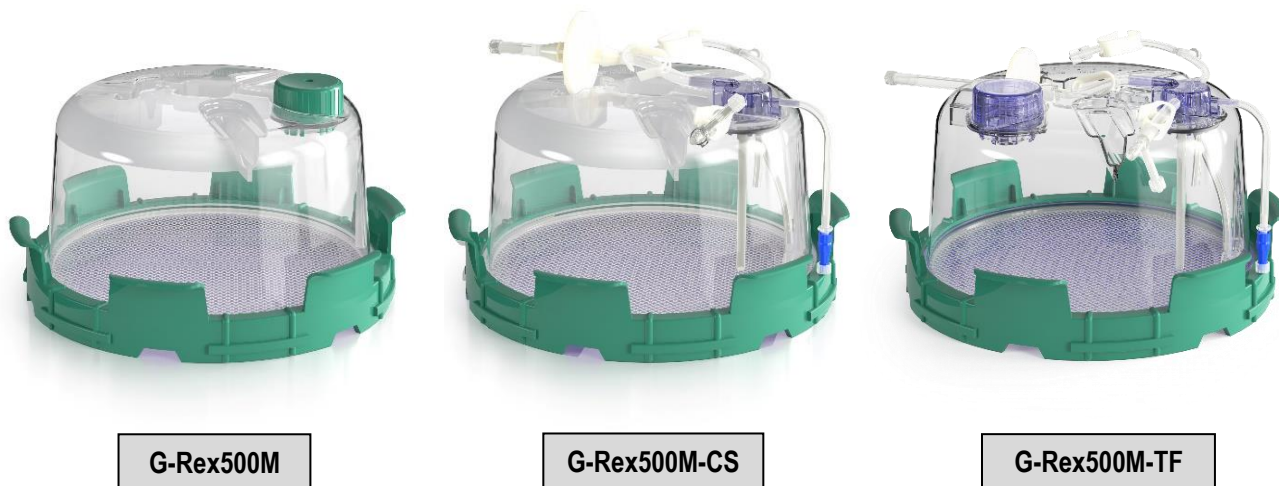

# WILSONWOLF

## TECHNICAL DATA SHEET



### 2) Key Characteristics & Dimensions

Specification	G-Rex500M	G-Rex500M-CS	G-Rex500M-TF
Functional Use	Open System	Closed System	
Cell Growth Area (membrane)	500 cm <sup>2</sup>		
Nominal Working Volume	5,000 mL		
Maximum Fluid Volume	5,500 mL		
Approximate Dry Weight	2.50 lbs (1,135 grams)	2.75 lbs (1,249 grams)	2.82 lbs (1,278 grams)
Estimated Filled Weight @ Nominal Working Volume	13.53 lbs (6,135 grams)	13.78 lbs (6,249 grams)	13.84 lbs (6,278 grams)
Maximum Product Height	7.13 inch (18.1 cm)	7.88 inch (20.0 cm)	7.88 inch (20.0 cm)
Maximum Product Diameter	12.27 inch (31.2 cm)		
Sample Tubing Line Connection	N/A (open system)	<ul style="list-style-type: none"> <li>▪ MicroClave® Connector</li> <li>▪ Total Sample Line Residual Volume (~1.5 mL)</li> </ul>	
Reduction Tubing Line Connection Options	N/A (open system)	<ul style="list-style-type: none"> <li>▪ Female MPC</li> <li>▪ Female Luer Lock</li> <li>▪ Weldable PVC Tubing (0.114" ID / 0.160" OD)</li> <li>▪ Weldable C-Flex Tubing (0.250" ID / 0.375" OD)</li> </ul>	
Harvest Tubing Line Connection Options	N/A (open system)	<ul style="list-style-type: none"> <li>▪ Male Luer Lock</li> <li>▪ Weldable PVC Tubing (0.114" ID / 0.160" OD)</li> </ul>	

# WILSONWOLF

## TECHNICAL DATA SHEET



### 3) Packaging

Feature	G-Rex500M	G-Rex500M-CS	G-Rex500M-TF
Packaging Materials	PETG Tray & Retainer, Tyvek® Lid	PETG Tray & Retainer, Tyvek® Lid	PETG Tray & Retainer, Tyvek® Lid
Devices per Package	1		
Minimum Shipping Quantity	2-pack (Case)		
Case Dimensions (Length x Width x Height)	L 18.38 in (46.7 cm) W 15.38 in (39.1 cm) H 15.63 in (39.7 cm)		
Case Weight	10.91 lbs (4.95 kg)	11.21 lbs (5.09 kg)	11.59 lbs (5.26 kg)
Shelf Life	3-years		
Recommended Storage Conditions	Ambient (64 - 82°F (18 to 28°C), No Humidity Controls)		


# WILSONWOLF

## TECHNICAL DATA SHEET




### 4) Packaging Images

#### Individual Packaging

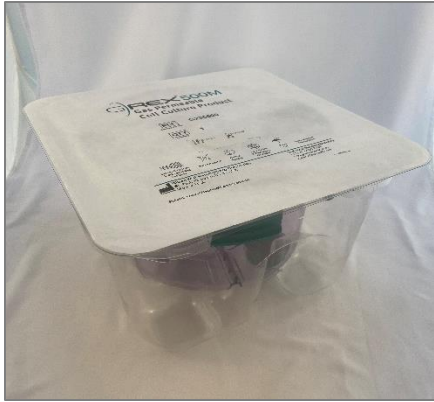

G-Rex500M

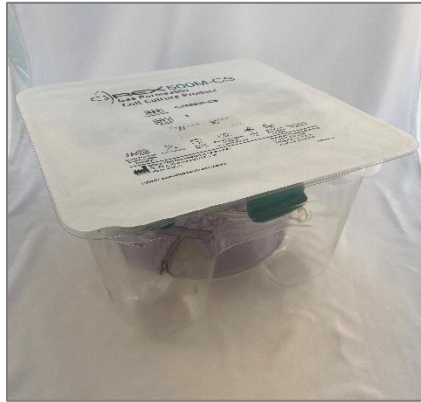

G-Rex500M-CS



G-Rex500M-TF

#### 2-Pack Cases



G-Rex500M



G-Rex500M-CS



G-Rex500M-TF

# WILSONWOLF

## TECHNICAL DATA SHEET



### 5) Other Product Information


Feature	Value
Certificate of Compliance (CoC)	Provided with Every Shipment
Viral Penetration Tested (Membrane Material only)	ASTM F1671
No Chemical Additives Statement	G-Rex Products Do Not Use Latex, BPA, Phthalates, Nitrosamines, Residual Solvents or Mold Release Agents in their Component Manufacturing or Final Device Assembly.
Electronic Instructions for Use (eIFU)	Visit <a href="http://www.wilsonwolf.com">www.wilsonwolf.com</a>

# WILSONWOLF

## TECHNICAL DATA SHEET



### 6) Product Materials

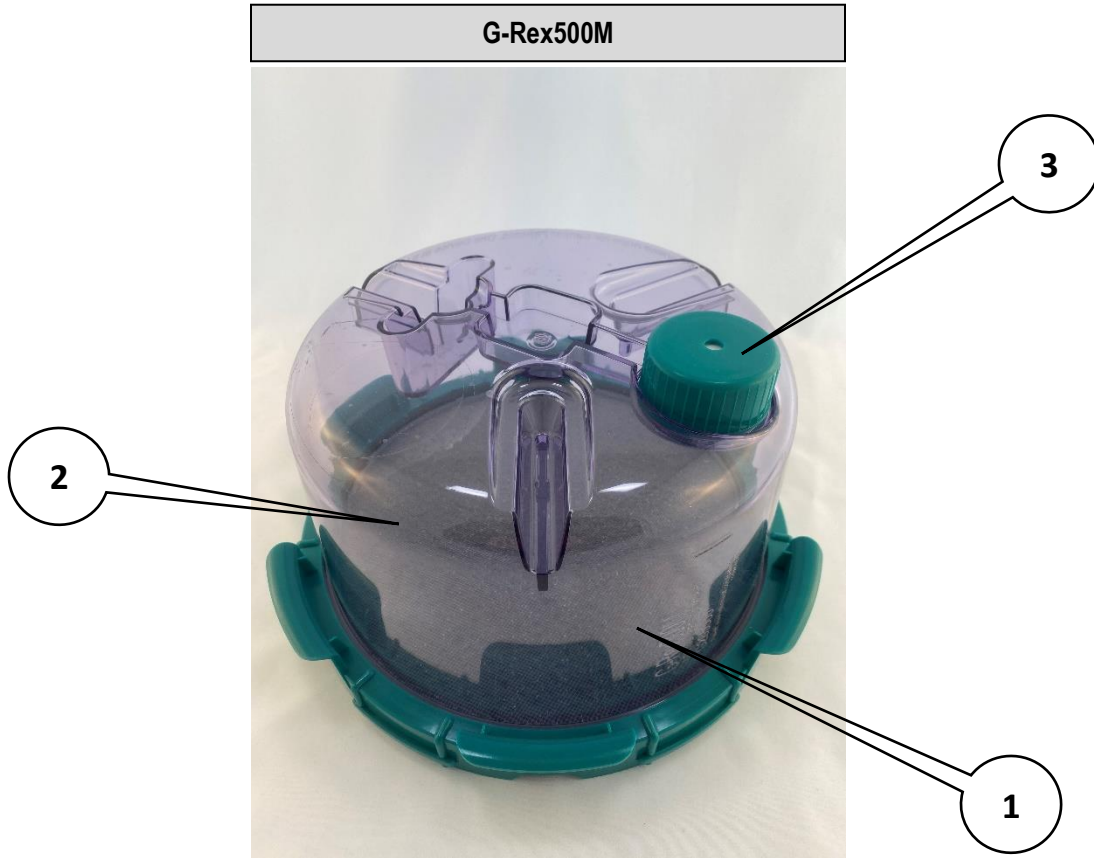


Fluid Pathway Components (numbered)

#	Component Name	Material	Material Grade
1	Gas Permeable Membrane	Silicone Elastomer	USP Class VI, ISO 10993
2	Shell (Vessel)	Polycarbonate	USP Class VI, ISO 10993
3	Cap with Vent Filter	Cap: Polypropylene, Filter Membrane: Hydrophobic ePTFE (0.2 μm)	USP Class VI, ISO 10993



# WILSONWOLF

## TECHNICAL DATA SHEET

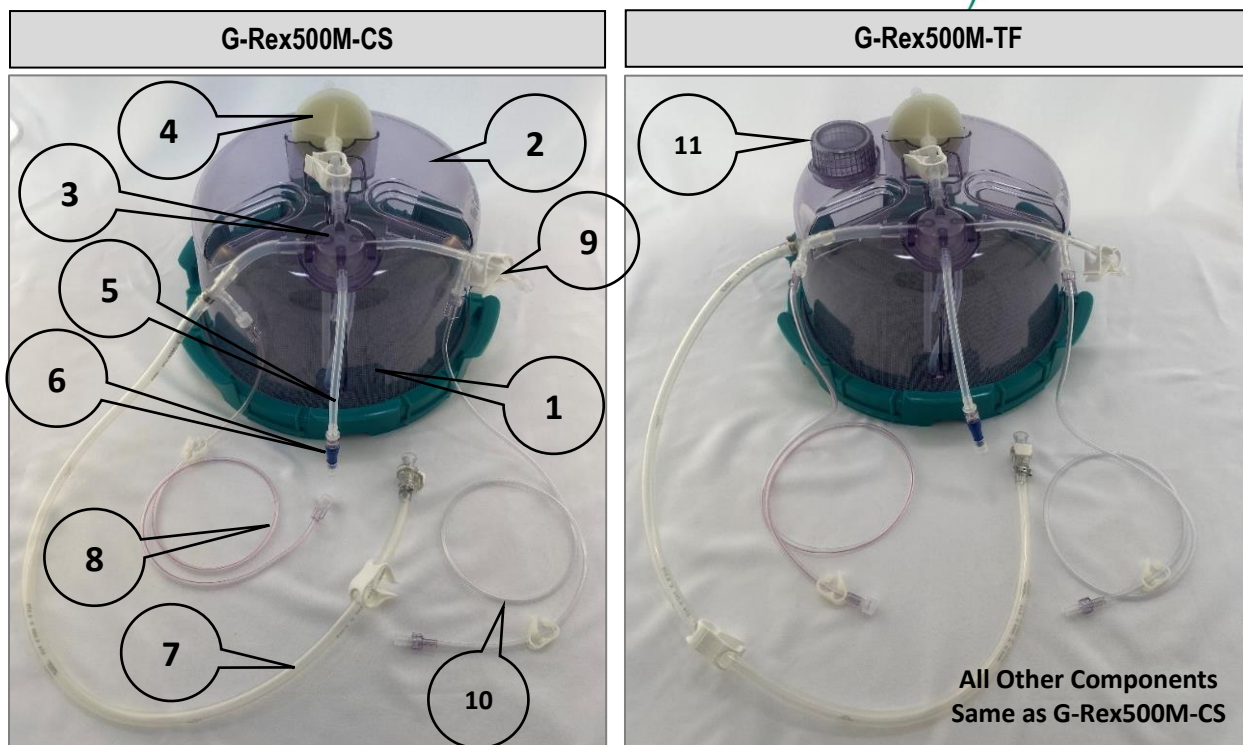

All Other Components Same as G-Rex500M-CS

Fluid Pathway Components (numbered)

#	Component	Material	Material Grade
1	Gas Permeable Membrane	Silicone Elastomer	USP Class VI, ISO 10993
2	Shell (Vessel)	Polycarbonate	USP Class VI, ISO 10993
3	Multi-Port Cap	Polypropylene	USP Class VI, ISO 10993
4	Pall® Emflon® II Vent Filter (non-fluid pathway)	0.2 µm, 50mm Diameter, Hydrophobic Polyvinylidene fluoride (PVDF), Emflon® II membrane in mini Kleenpack™ filter capsule (Pall® Product No. KM5V002P2G100)	USP Biological Reactivity Tests, In Vivo <88>
5	Sample Line Tubing	Silicone Elastomer	USP Class VI, ISO 10993
6	MicroClave® Connector	Body: Polybutylene Terephthalate (PBT), Seal: Silicone, Spike: Polycarbonate, Lubricant: Silicone.	FDA Class 2 Device USP Class VI, ISO 10993
7	Reduction Line Tubing	C-Flex®, Thermoplastic Elastomer	USP Class VI, ISO 10993
8	Weldable Reduction Line	Polyvinyl chloride (PVC), non-toxic, non-phthalate	USP Class VI, ISO 10993
9	Harvest Line Tubing	Silicone Elastomer	USP Class VI, ISO 10993
10	Weldable Harvest Line	Polyvinyl chloride (PVC), non-toxic, non-phthalate	USP Class VI, ISO 10993
11	Cap (non-vented)	Polycarbonate	USP Class VI, ISO 10993
-	Various Luer Fittings	Polycarbonate	USP Class VI, ISO 10993
-	Various Luer Caps	Polypropylene	USP Class VI, ISO 10993

Wilson Wolf Manufacturing LLC

St. Paul, MN 55112 USA

(651) 628-9259

info@wilsonwolf.com

www.wilsonwolf.com

MM 0006 Rev B

**7) Regulatory / Compliance**

Company Certifications	ISO 13485
OUS Product Classification	GMP raw material for the manufacturing of cellular products, also known as an ancillary material per ISO/TS 20399-3:2018(en) Biotechnology - Ancillary Materials Present During The Production Of Cellular Therapeutic Products - Part 3: Best Practice Guidance For Ancillary Material Users, and USP <1043> Ancillary Materials for Cell, Gene, and Tissue-Engineered Products.
FDA Est. Reg. # and Est. Type	3004982543, Medical Device Manufacturer
FDA Product Classification	Device Class: I Classification Name: System, Suspension, Cell Culture Product Code: KJF Regulation Number and Description: 21 CFR 864.2240 Cell and Tissue Culture Supplies and Equipment Regulation Medical Specialty: Hematology Premarket Submission Type: 510(k) Exempt
Patents	<a href="http://www.wilsonwolf.com/patents">www.wilsonwolf.com/patents</a>

**8) Legal Disclaimer**

Wilson Wolf Manufacturing LLC makes no claims regarding the performance of this product for clinical treatment or therapeutic applications. It is the responsibility of the end user to assess its suitability for specific applications.